

## ■ WARNING INDICATIONS

### OVERLOAD INDICATION

Scale overload. (max. weight capacity is indicated on the scale).  
Please step off to avoid damage.




### LOW BATTERY INDICATION

Battery power is low. Please install new battery.



### RE-TEST

Display "[-]" when you step off the scale, means this weight reading with error. Please re-test.

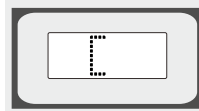

## ⚠ ADVISE FOR USE & CARE

- 1.The glass platform will be slippery when it is wet. Keep it dry!  
Keep the scale in cool place with dry air.
- 2.Keep standing still during the course of weighing. DO NOT strike, shock or drop the scale.
- 3.Treat the scale with care as it is a precise instrument. Always keep your scale in horizontal position.
- 4.Clean the scale with a damp cloth and avoid water getting inside the scale. DO NOT use chemical/abrasive cleaning agents.
- 5.If failed to turn on the scale, check whether battery is installed, or battery power is low. Insert or replace new battery.
- 6.If there is error on display or can not turn off the scale for a long time, please take off the battery for about 3 seconds, then install it again to remove soft-fault. If you can not resolve the problem, please contact your dealer for consultation.
- 7.Do not use it for commercial purpose.

IMPORTED & MARKETING BY:

**Passim Medichem Agencies**

#186, Ind. Area, Chandigarh 160 002

Corporate Office:

#902, 9th Floor, WallFort House, S.V.Road

Near Citi Mall, Gaegaon (W)


Mumbai 400 062, INDIA

Customer care:

Website: [www.drodin.in](http://www.drodin.in)

Email ID: [customercare@drodin.in](mailto:customercare@drodin.in)

Contact: +91 9878785333, 9359490504

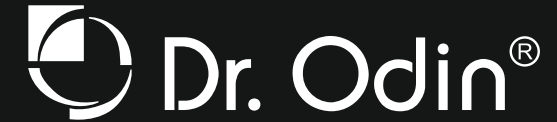

# Electronic Personal Scale


## MODEL NO.: EB-7010

## Instruction for use



ISO 9001 CERTIFIED BY SGS  
EB-36 SERIES

## ■ TO REPLACE BATTERY

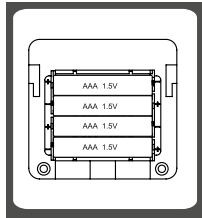


4x1.5V (AAA)batteries.

If you want to replace the batteries,

- A. Open the battery cover at the bottom of scale;
- B. Take out the used batteries;
- C. Install new battery by putting one side of the battery, and then pressing down the other side. Please take note to the polarity.

## ■ LED DISPLAY

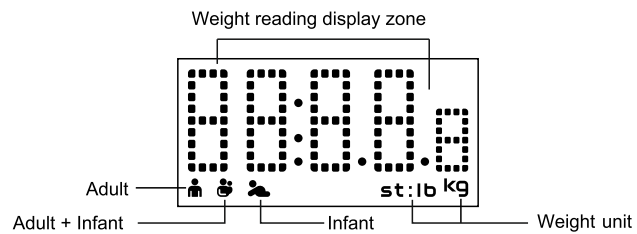

Remind that Adult should step on the scale, and what the LED shows is the weight for Adult.




Remind that Adult should step on the scale with infant.



Remind that what the LED shows is the weight for infant.

Remark: The above three icons will be only available under the Infant Mode.

## ■ BUTTON FUNCTION

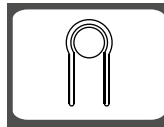

Button at the bottom of scale

When the scale is turned on:

1. Press for a white: Convert the mode between Adult Mode and Infant Mode
2. Press shortly : Convert the weight unit among kg / lb /st


Remark: The unit system of scale depends on sales destination or the requirement of distributor.

## ■ OPERATION

### Adult Mode: For adult weighing



Step 1: To make sure accuracy, place the scale on a hard, flat floor and avoid carpet or soft surface.

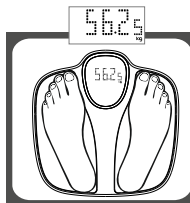

Step 2: Adult steps on the scale gently, the scale will be on automatically. Stand for a while and wait until the LED shows the weight reading. The reading will flash three times and locked (as shown like the picture).


### Infant Mode: Weighing adult & Infant respectively




Step 1: To make sure accuracy, place the scale on a hard, flat floor and avoid carpet or soft surface.



Step 2: Adult steps on the scale gently, the scale will be on automatically. Stand for a while and wait until the LED shows the weight reading. Once the reading is stable, the LED will show Zero, then the "Adult+Infant" icon will flash (as shown like the picture)

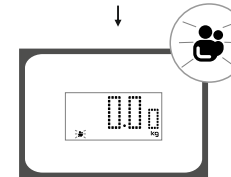

Step 3: Adult steps off the scale and steps on again with infant. Stand for a while and wait until the LED shows the weight of infant. The LED will show the weight of adult and infant respectively for three times.

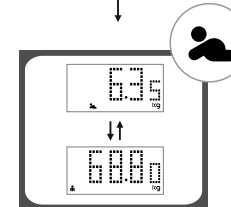

Remark: the scale will be off automatically within seconds if there is no operation.