

DISCLAIMER OF WARRANTY

This warranty shall not apply if this product is damaged due to accident, misuse, neglect or other causes not related to any manufacturing defects. Any applicable implied warranties including warranties of merchantability and fitness for a particular purpose and hereby expressly disclaimed.

Proof of purchase in the form of a bill or sale or receipted invoice which is evidence that the product is within the warranty period must be presented to obtain warranty services.

NOTE

THE CUSTOMER WILL HAVE TO SHIP THE PRODUCT TO THE COMPANY AT IT'S OWN EXPENSE TO AVAIL WARRANTY SERVICES.

WARRANTY CARD	
Model No.....	Lot No.....
Invoice No.....	Date of Purchase.....
Purchased by :.....	Contact No. :.....
Address.....	
Dealer Name.....	
Dealer Sign & Stamp	

Manufactured by:

Passim Lifesciences Ltd.

Plot No. 45, Ind. Area, Phase-II,
Panchkula - 134113 (Haryana) - **INDIA**

Test Lic. No.: MFG/TL/MD/2023/000112

Marketed by:

Odin Healthcare Pvt. Ltd.

Plot No. 45, Ind. Area, Phase-II,
Panchkula - 134113 (Haryana) - **INDIA**

**For any Complaint/ Suggestion
please contact :**


Customer Care No.: 1800 309 3009

(Timing: 9am -7pm, Mon. - Sat.)

Dr. Odin[®]

Single Bed Warmer

OBW101

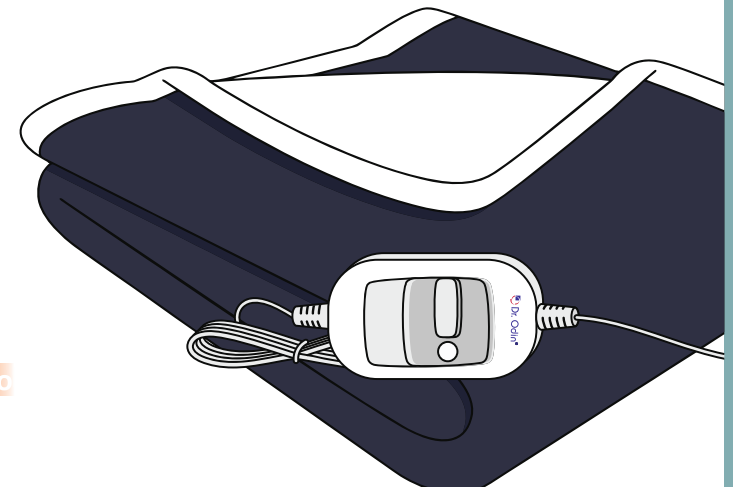

- ❖ 70 V
- ❖ 3 Heat
- ❖ Dual Safety
- ❖ Overheat Protection
- ❖ 12 Hours Auto

IMPORTANT SAFEGUARDS

In order to reduce the risk of fire, electrical shock or injury to person while using the electric appliance, it is important to follow these basic safety precautions:

- Intend to use ONLY as a bed warmer. This is not an over blanket.
- NOT intended for use with pets.
- Unfold completely for proper use.
- Check bed warmer frequently to ensure it is NOT bunched or folded as this could cause excel heating.
- If the bed warmer does not function properly, discontinue use immediately and contact our customer care.
- Regularly inspect the cord, plug, power supply and bed warmer for damages. If damaged, discontinue use immediately and contact customer care
- Bed warmer should not hangover the edges of the bed.
- Do not leave disabled persons, infants or young children unattended while in use.
- Long applications on the high heat setting can result in skin irritation or burns.
- Keep controls unit away from open windows.
- Do not pierce with sharp objects.
- Use bed warmer only on an alternate current (AC) supply circuit of correct voltage.
- Do not use the bed warmer on the power supply generated by generators.

CARE & MAINTENANCE

• DO NOT DRY- CLEAN OR MACHINE WASH

- Cleaning liquids can deteriorate the insulation and covering of the heating unit.

• DO NOT IRON

• SPOT CLEANING

Sponge the soiled are using luke warm water and mild detergent or soap. Follow with a sponging of clean water to remove excess water, compress the wet area between absorbent towels.


INSTRUCTION FOR USE

- Spread the bed warmer under the bed sheet and over the bed mattress. No portion of the bed warmer should be put under the pillow or tucked under the mattress.
- Make sure the bed warmer is spread smoothly without any folds.
- The control unit of the bed warmer should be at the head side of the bed and should be with in the easy reach to adjust the heat settings.
- After spreading the bed warmer, cover the bed with a blanket or a quilt for enhanced

warmth, so that the heat which is being generated gets accumulated inside the blanket/ quilt and no heat is lost in open air.

- Connect the bed warmer to the main power.
- For initial warming of the cold bed, set the heat setting to level “3” (Red indication LED) in order to provide full heating.
- Allow 20-25 minutes for pre-heating of the bed to achieve luxurious warmth.
- Usage of bed warmer on generators or on high voltage can damage the bed warmer.

OPERATION & CONTROLS

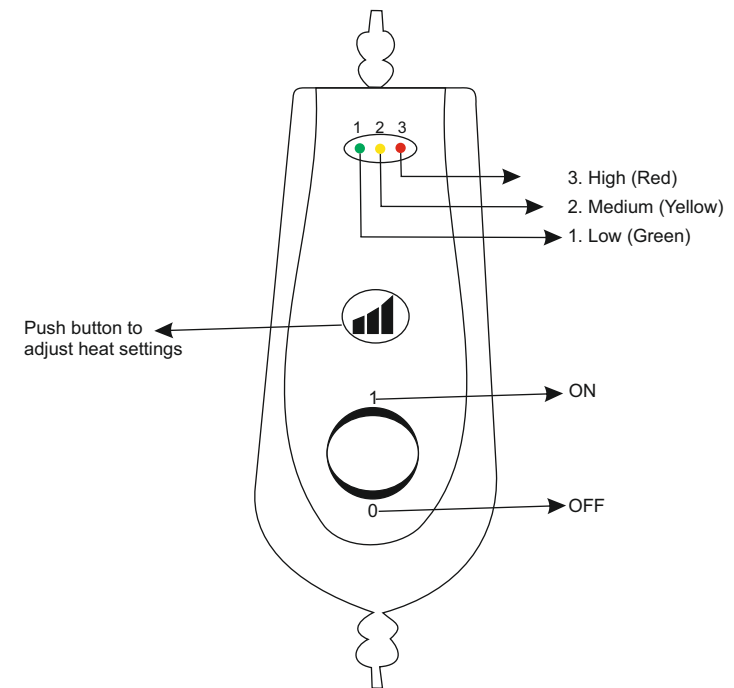

WARRANTY: LIMITED ONE YEAR WARRANTY

Dr. Odin Electric Bed Warmer warrants to the original purchaser that this product shall be free from defects in material, design and craftsmanship for a period of one year from the date of purchase (The Warranty Period). If this product is determined to be defective during the warranty period, we agree to either repair or replace, at its option. Contact customer care to obtain our step to step instructions to receive, repair/ replacment services.